


BIOREACTOR TUBING KITS

SINGLE-USE & CUSTOMIZABLE FOR BENCHTOP BIOREACTORS

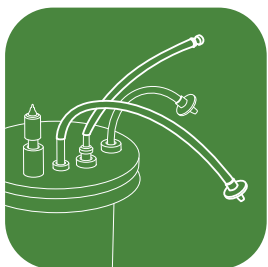

SUS
www.advantapure.com

Bioreactor Tubing Kits are built-to-order and clean room packaged. The Single-Use kits are custom manufactured and assembled from a variety of components — tubing, filters, connectors — to offer faster turnaround time and increased processing capacity. Made from platinum-cured silicone or AdvantaFlex® biopharmaceutical grade TPE tubing, AdvantaPure's Bioreactor Tubing Kits provide convenience and greater efficiency.

KEY FEATURES

- Components meet relevant ISO and USP standards, including Class VI
- Assembled and packaged in ISO 7 (Class 10,000) clean rooms
- AdvantaFlex TPE tubing is sterile weldable and heat sealable
- Tubing is certified free of animal-derived ingredients
- Documented quality control; lot traceable with identification on bags
- Tubing material validation and extractables test portfolios available upon request
- Optional RFID tags easily attach to assemblies for identification purposes
- Compounds used to make tubing have a Master File with the U.S Food and Drug Administration (APST Silicone Drug Master File #26598; AdvantaFlex TPE Drug Master File #28810)

BIOREACTOR TUBING KITS

SINGLE-USE & CUSTOMIZABLE FOR BENCHTOP BIOREACTORS



AdvantaPure's Bioreactor Tubing Kits:



- Offer faster turnaround time for your bioreactor and increase your processing capacity
- Decrease cross contamination risks
- Reduce repetitive motion injuries from assembly processes
- Minimize the number of inventoried components you need to stock
- Provide a single, animal-derived-ingredient-free material contact surface (both tubing and molded components)
- Design your ready-to-use assembly to your specifications


SUS

■ Kits with AdvantaFlex® TPE tubing allow for easy welding and sealing.



Bioreactor Tubing Kits are Custom Assembled for Your Process

- Specify your tubing material and length for addition, vent, sparge and harvest lines.

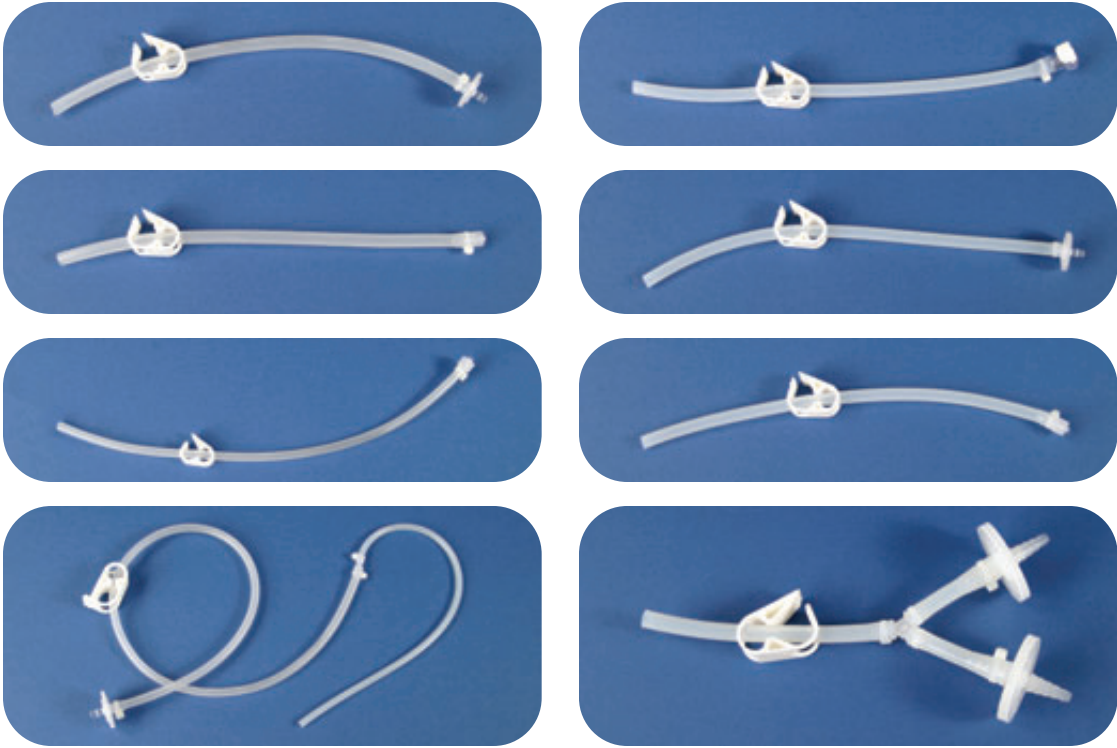

- You identify the tubing sizes, materials and lengths, plus filters and connectors.

SUS



- AdvantaFlex TPE lines are available with sealed ends; silicone lines are available plugged.



BIOREACTOR TUBING KITS

SINGLE-USE & CUSTOMIZABLE FOR BENCHTOP BIOREACTORS



Select from AdvantaFlex TPE Tubing or APST Silicone Tubing



 **AdvantaFlex**[®]
BioPharmaceutical Grade Tubing



AdvantaFlex addresses the need for a flexible, sterilizable, moldable, weldable and heat sealable biopharmaceutical tubing without silicone oils. AdvantaFlex maintains its physical properties following sterilization processes, remains translucent for visible product flow and does not become gummy. The tubing resists kinking and is an excellent component for Single-Use systems.

Learn more at www.advantapure.com/advantaflex.htm



APST UNREINFORCED SILICONE TUBING

Low-volatile grade, platinum-cured silicone tubing is cleanroom produced for critical pharmaceutical applications. APST has undergone extensive physical, chemical, and biological testing and meets USP Class VI, FDA CFR 177.2600, ISO 10993, European Pharmacopoeia 3.1.9, and 3-A standards.



Learn more at www.advantapure.com/apst-silicone-tubing.htm

SUS

Address: NewAge Industries AdvantaPure[®]
145 James Way
Southampton, PA 18966 U.S.A.

 **AdvantaPure**[®]
NEWAGE INDUSTRIES



DISTRIBUTED BY:

Toll Free Phone: **888-755-4370**
Phone: **215-526-2151**
Toll Free Fax: **888-258-4293**
Fax: **215-526-2167**
Web Site: **www.advantapure.com**
E-mail: **sales@advantapure.com**

©2016 NewAge[®] Industries, Inc. Specifications are subject to change without notice. AdvantaPure products are produced by NewAge[®] Industries, Inc. AdvantaFlex[®], Purity in Fluid Flow Systems[®], NewAge[®], NewAge Industries AdvantaPure[®], and the "N" logo are registered trademarks of NewAge[®] Industries, Inc.



ESOP Employee Owned
for Your Benefit

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term NewAge.