

DocuLink® is a service designed to save time and effort during document retrieval. DocuLink identifies a single use system on the product itself—not only on the bag from which it was removed—by use of a permanently-affixed label encapsulated into the molded silicone connector.

DocuLink lists the individual part number and lot number to allow quick lookup of documents such as certificates of compliance and certificates of analysis.




KEY FEATURES

- Visually identifies each manifold with a part number and/or lot number
- Label is permanently molded into the silicone connection
- Does not come into contact with the critical stream
- Allows easy tracking of a manifold assembly back to its origins
- Handles clean-in-place (CIP), steam-in-place (SIP), gamma irradiation, and autoclave processes
- Cannot fall off like adhesive labels
- Molded and packaged in an ISO Class 7 (Class 10,000) clean room
- Identification can be linked to an online data management site
- Documented quality control; lot traceable with identification on polybags




IDENTIFICATION
www.advantapure.com

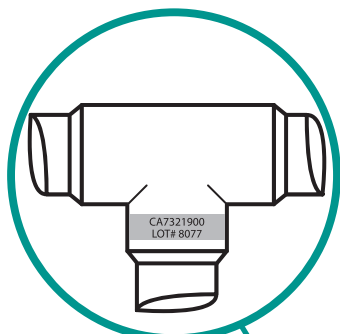


Customize to fit your needs

- Every individual manifold is manufactured with a permanently affixed label including part number and lot number
- Enables tracking of manifold assembly back to its origins
- Permanently molded into the silicone connector
- Molded and packaged in a Class 7, (Class 10,000) ISO-certified clean room
- Handles clean-in-place (CIP) or steam-in-place (SIP) processes
- Autoclavable
- Gamma sterilizable

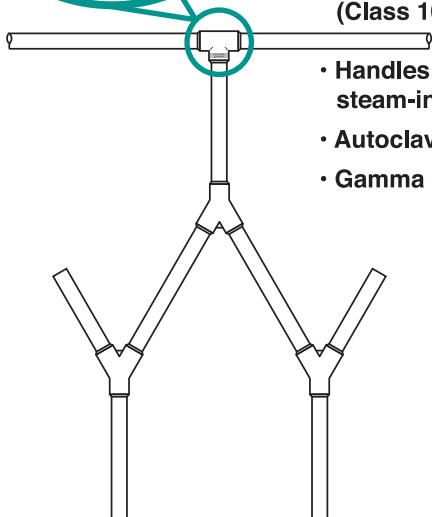

Manifold I.D. Number can be linked to an online data management site (optional)

- Fast, easy access to certifications for each individual manifold & its AdvantaPure-assembled components
- Facilitates gathering of documents for FDA audits
- A less expensive method of document storage & retrieval

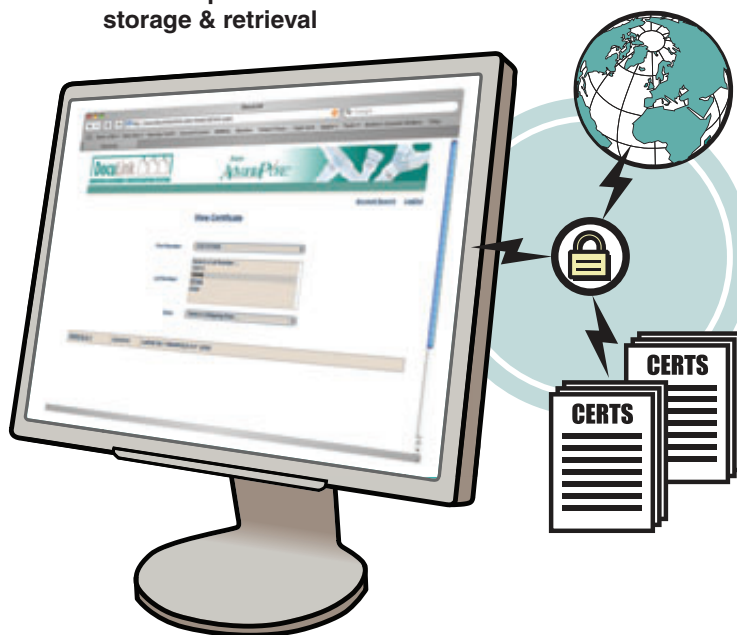

IDENTIFICATION

IDENTIFICATION



• Login to secure Web site



• Select manifold part number from drop down box




• Select lot number from list to display available documents




• View or download PDF documents



OTHER IDENTIFICATION METHODS

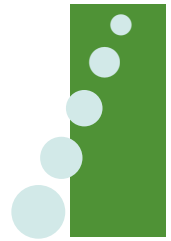

GammaTAG® GAMMA STERILIZABLE RFID TAGS

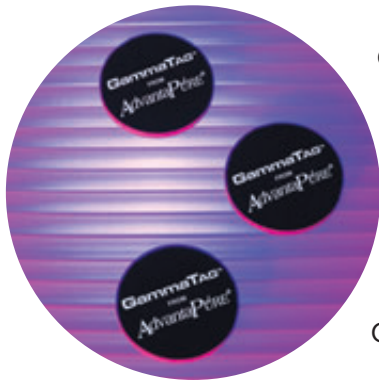

GammaTag is the first read/write RFID tag that's safely sterilizable by gamma radiation up to 45 kGy. Use GammaTag to reliably identify, track, and store electronic data of Single-Use bioprocess components from inception to disposal.

GammaTag is available from NewAge® Industries' Verigenics® group. Visit www.verigenics.com.

Gamma irradiation services of complete systems, with or without GammaTags, are also available.

GammaTag was developed in partnership with **ProcessHQ**
PURE AND SIMPLE™

Advanta LABEL® HOSE LABELING SYSTEM



Clear identification, cleanliness, sterilizability, and permanence are offered with AdvantaPure's AdvantaLABEL® tagging system.

AdvantaLABEL eliminates the issues involved with old-style metal, plastic, or even paper tags. Since the label is sealed to the hose in clean silicone, it can be either steam sterilized (SIP) or chemically cleaned (CIP) many times without degradation of the label.

Use AdvantaLABEL to denote any information you deem critical – part number, date of manufacture, planned replacement date, more. Contact your AdvantaPure sales representative or visit www.advantapure.com/advantalabel.htm for information.



IDENTIFICATION

Purity in Fluid Flow Systems®




Address: NewAge Industries AdvantaPure®
145 James Way
Southampton, PA 18966 U.S.A.



Toll Free Phone: **888-755-4370**
Phone: **215-526-2151**
Toll Free Fax: **888-258-4293**
Fax: **215-526-2167**
Web Site: **www.advantapure.com**
E-mail: **sales@advantapure.com**

DISTRIBUTED BY:

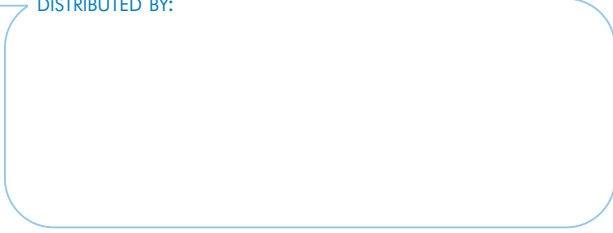

©2016 NewAge® Industries, Inc. Specifications are subject to change without notice. AdvantaPure products are produced by NewAge® Industries, Inc. AdvantaLabel®, DocuLink®, GammaTag®, Purity in Fluid Flow Systems®, Verigenics®, NewAge Industries AdvantaPure®, NewAge®, and the "N" logo are trademarks of NewAge® Industries, Inc.

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term New-Age.

