


OVERMOLDED TRI-CLAMP® ENDS

OF HYTREL® & POLYPROPYLENE

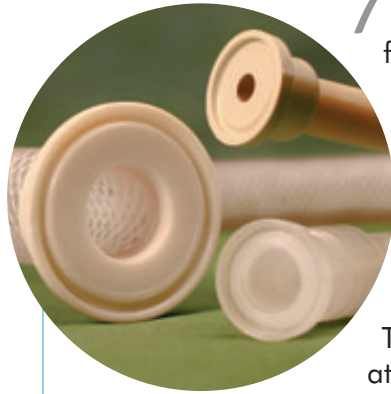

An alternative for sanitary connections is now available from AdvantaPure . . . Overmolded Tri-Clamp Ends for PVC, TPE, and TPR tubing and hose.

Injection molded from USP Class VI Hytrel® or polypropylene, Overmolded Tri-Clamp Ends feature Tri-Clamp and mini Tri-Clamp styles up to 1-1/2".

Overmolded Tri-Clamp Ends greatly reduce bioburden entrapment situations associated with barb-type fittings. They cost less than silicone counterparts, making them an attractive option for pharmaceutical, medical, surgical, food, beverage, and dairy applications.




KEY FEATURES

- Injection molded from USP Class VI Hytrel and polypropylene materials
- Certified free of animal-derived ingredients
- Ends are permanently attached using a thermal process — no leachable adhesives are involved
- Smooth interior for excellent flow and performance
- Reduces entrapment issues
- Provides a sanitary connection for thermoplastic tubing and hose
- Excellent component for single-use applications
- Less costly than similar components of silicone
- Available in both Tri-Clamp and mini Tri-Clamp
- Sizes range from mini through 1-1/2" Tri-Clamp
- Easily sterilized/cleaned
- Hytrel molds directly to PVC; polypropylene molds directly to TPR and TPE tubing, such as AdvantaPure's AdvantaFlex®
- Assemblies made to order
- Documented quality control




HYTREL®/PP
www.advantapure.com

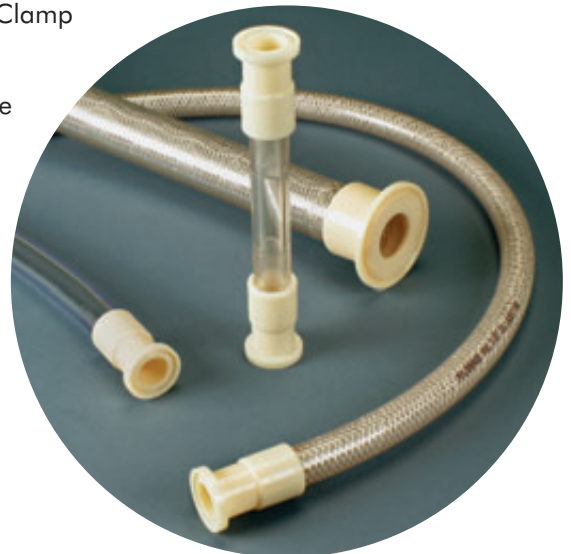

OVERMOLDED TRI-CLAMP® ENDS

OF HYTREL® & POLYPROPYLENE

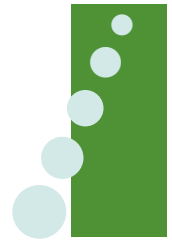

SIZES

Tubing/Hose I.D. (Inch)	Tubing/Hose O.D. (inch)	End Style
1/8	1/4	Mini Tri-Clamp
1/4	3/8	Mini Tri-Clamp
1/4	7/16	Mini Tri-Clamp
1/4	1/2	Mini Tri-Clamp
3/8	5/8	Mini Tri-Clamp
1/2	3/4	Mini Tri-Clamp
5/8	7/8	Mini Tri-Clamp
1/4	3/8	1" or 1-1/2" Tri-Clamp
1/4	7/16	1" or 1-1/2" Tri-Clamp
3/4	1-1/8	1" or 1-1/2" Tri-Clamp
1	1-3/8	1-1/2" Tri-Clamp



HYTREL®/PP

Purity in Fluid Flow Systems®



Made in U.S.A.

Address: NewAge Industries AdvantaPure®
145 James Way
Southampton, PA 18966 U.S.A.



DISTRIBUTED BY:

Toll Free Phone: **888-755-4370**
Phone: **215-526-2151**
Toll Free Fax: **888-258-4293**
Fax: **215-526-2167**
Web Site: **www.advantapure.com**
E-mail: **sales@advantapure.com**

©2015 NewAge® Industries, Inc. Specifications are subject to change without notice. AdvantaPure products are produced by NewAge® Industries, Inc. Purity in Fluid Flow Systems®, NewAge Industries AdvantaPure®, NewAge®, and the "N" logo are trademarks of NewAge® Industries, Inc. Hytrel® reg TM DuPont. • Tri-Clamp® registered trademark Alfa Laval Inc.



ESOP Employee Owned
for Your Benefit

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term New-Age.