


# APEWF WIRE REINFORCED EPDM HOSE

WHITE EPDM CORE; HIGHLY FLEXIBLE

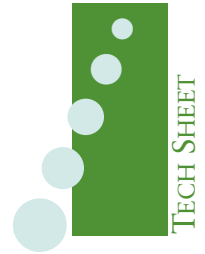

Premium-quality hose for elevated pressure and temperature applications in the pharmaceutical, food, and beverage processing industries is available from AdvantaPure. APEWF is specially engineered to handle a wide range of fluids and chemicals.

APEWF is reinforced with a wire helix, making it suitable for suction or discharge applications. The hose offers flexibility and a smooth, sanitary inner surface for unrestricted flow for uses such as washdown stations, coolant transfer, glycol lines, and others.



## KEY FEATURES

- Excellent chemical resistance properties make this hose ideal for many demanding applications
- Designed for suction or discharge applications
- Helical wire reinforcement helps withstand a higher level of internal pressures and maintain stated vacuum ratings
- Sanitary, smooth, white EPDM rubber core resists bacteria build-up and complies with FDA standards for use in food and beverage applications such as wine and beer, potable water, and fruit juices
- Odorless and tasteless core – fluids remain pure
- Handles a more extreme range of working temperatures than standard EPDM hoses
- Offers improved flexibility over similar fluoropolymer hoses
- Ideal for handling CIP solutions; may also be steam sterilized up to 302°F (150°C) for a maximum of 30 minutes
- Custom designed hose construction for flexibility – easily handles bends, turns, corners
- Abrasion-resistant smooth rubber cover
- Sizes range from 3/4 in. (19 mm) up to 3 in. I.D. (76 mm)
- Full line of crimp fittings and custom assemblies offered
- Add stainless steel armor casing or cuffs for an extra layer of protection from crushing, kinking, and harsh environments while still allowing flexibility
- Laser-etched stainless steel collars for quick, visual hose identification available
- Hose Track® RF Identification and Lifecycle Analysis System available – RFID tags uniquely identify individual hose assemblies, create audit trails, and track key wear-related events


EPDM RUBBER  
[www.advantapure.com](http://www.advantapure.com)



PURITY IN FLUID FLOW SYSTEMS®

# APEWF WIRE REINFORCED EPDM HOSE

WHITE EPDM CORE; HIGHLY FLEXIBLE

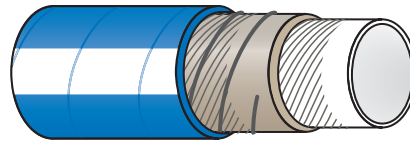

## SPECIFICATIONS

Product Number	Nominal I.D.		Nominal Wall		Nominal O.D.		Vacuum Rating		Working Pressure		Bend Radius		Weight (lbs./ft.)
	(in.)	(mm)	(in.)	(mm)	(in.)	(mm)	(in. Hg)	(Bar)	at 68°F (20°C)	(PSI)	(Bar)	(in.)	
APEWF-W-0750*	.75	19	.236	6	1.22	31	26.58	.9	450	32	2	66	.46
APEWF-W-1000*	1	25	.236	6	1.46	37	26.58	.9	400	28	3	88	.56
APEWF-W-1500*	1.5	38	.256	6.5	2.01	51	26.58	.9	300	20	5	135	.94
APEWF-W-2000*	2	51	.275	7	2.56	65	26.58	.9	225	16	7	200	1.28
APEWF-W-3000*	3	76	.295	8	3.58	91	23.62	.8	175	12	12	320	2.08

\*Limited stock availability – lead times and minimums may apply; call for details.

### NOTES:

Variables that impact working pressure include the type of fluid being conveyed through the hose, its temperature, pressure, and time. Field testing is recommended. Burst pressure is greater than or equal to three times the working pressure.

Temperature range: -40°F (-40°C) up to 248°F (120°C), depending on the product being conveyed. For sterilization purposes: up to 302°F (150°C) maximum for 30 minutes.

When products are used as part of an assembly, the pressure ratings of fittings may be less than hose pressure ratings above. Please consult your sales representative.

EPDM RUBBER

## Purity in Fluid Flow Systems®

Address: NewAge Industries AdvantaPure®  
145 James Way  
Southampton, PA 18966 U.S.A.



DISTRIBUTED BY:

Toll Free Phone: **888-755-4370**  
Phone: **215-526-2151**  
Toll Free Fax: **888-258-4293**  
Fax: **215-526-2167**  
Web Site: **www.advantapure.com**  
E-mail: **sales@advantapure.com**

©2020 NewAge® Industries, Inc. Specifications are subject to change without notice. AdvantaPure products are produced by NewAge® Industries, Inc. Hose Track®, Purity in Fluid Flow Systems®, NewAge Industries AdvantaPure®, NewAge®, and the "N" logo are trademarks of NewAge® Industries, Inc.



**ESOP** Employee Owned  
for Your Benefit

WARRANTY: NewAge Industries AdvantaPure products are warranted to be free from defects in material and workmanship. Any product found to carry such defects will be replaced without charge. This warranty is limited only to the replacement of NewAge Industries AdvantaPure products, item for item, which after examination by NewAge Industries AdvantaPure are deemed defective. NewAge Industries AdvantaPure specifically disclaims any other liability. The data supplied is provided as a helpful guideline and is believed to be reliable; however, nothing stated shall constitute a guarantee, recommendation, or warranty for any application.

NOTICE: NewAge Industries, Inc. has been a privately owned corporation since 1954. Neither the corporation nor its principals are affiliated, in any way, with any other group or entity either representing themselves as or being referred to by the term New-Age.